
LiNA Versascope™ Outflow Cannula

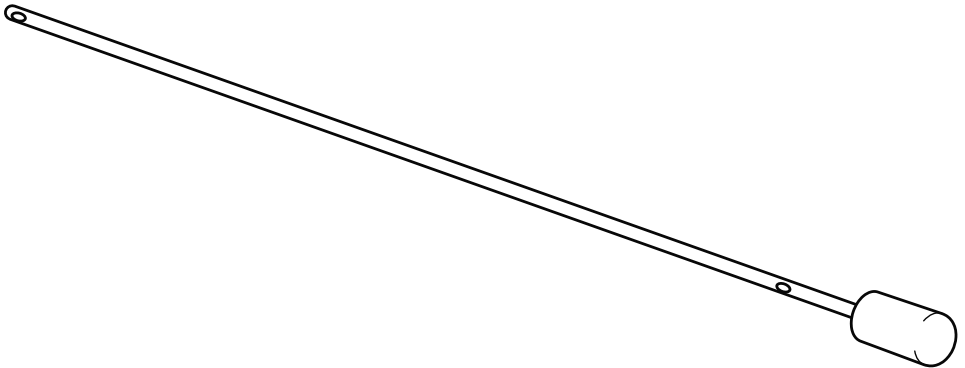
STERILE EO

CE
2274

LiNA Versascope™ Outflow Cannula

Ref: GMS720

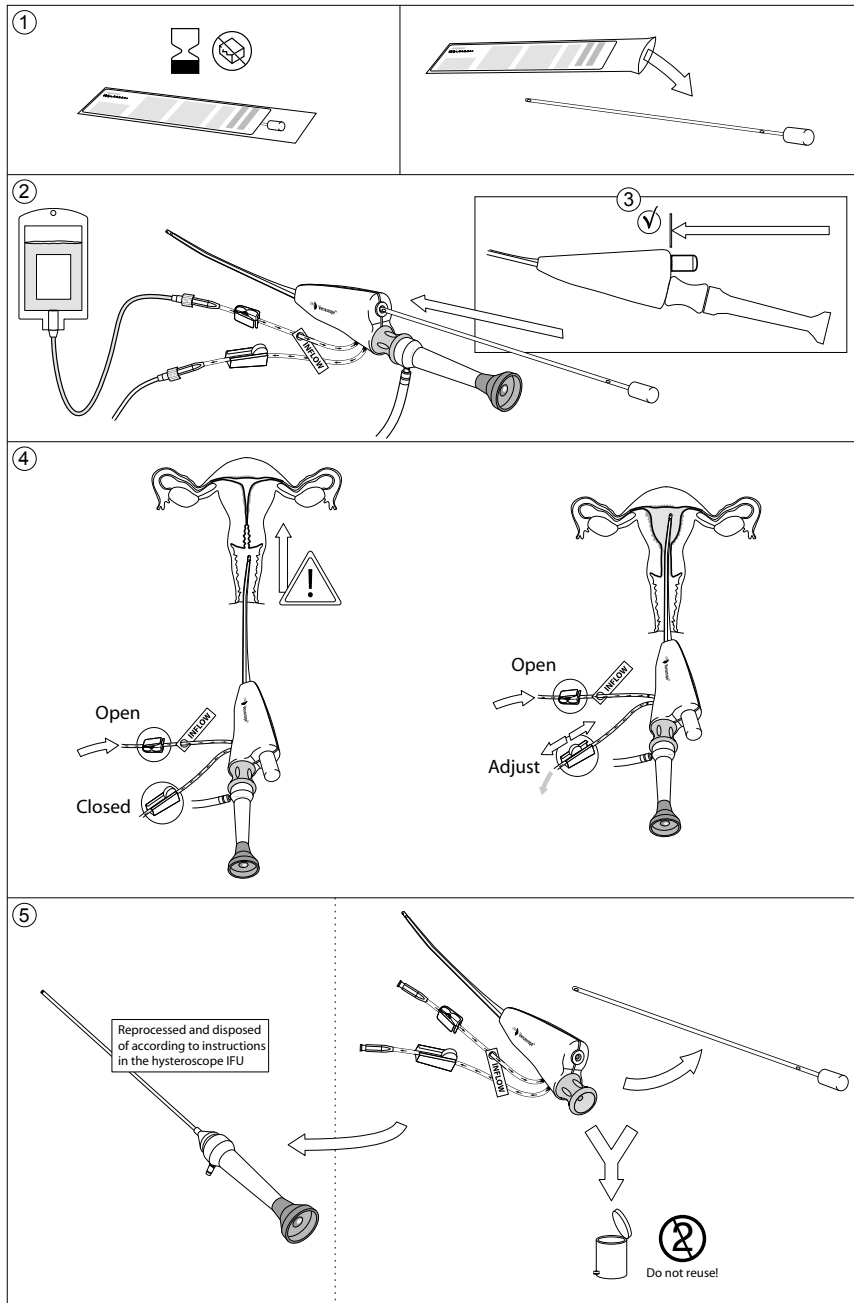Device description

LiNA Versascope™ Outflow Cannula GMS720 is a flexible, small diameter cannula intended to be used only with the LiNA Versascope™ Sheath. The LiNA Versascope™ Outflow Cannula is a single use, sterile device comprised of plastic components. The LiNA Versascope™ Outflow Cannula opens the expandable channel of the LiNA Versascope™ Sheath to provide for continuous flow, where fluids media is utilized for distension.

Intended use

The LiNA Versascope™ Outflow Cannula is used to open the expandable channel of the LiNA Versascope™ Sheath to provide the continuous flow, where the fluid media are utilized for uterine cavity distension.

Patient population

Adult females undergoing diagnostic or therapeutic hysteroscopy with the Versascope Schaft.

Contraindications

This product is not intended for use when hysteroscopic procedures are contraindicated.

Potential complications

When using a fluid distention medium, strict fluid intake and output monitoring is required. Excessive intravasation of distention fluid can lead to fluid overload. Potential fluid-related complications of Continuous Flow Hysteroscopy are:

- Hyponatremia,
- Hypothermia,
- Pulmonary edema,
- Heart Failure,
- Electrolyte Imbalances,
- Cerebral edema,
- Encephalopathy.


Other potential complications of Hysteroscopy include:

- Injury of the uterine wall/uterine perforation resulting in possible injury to bowel, bladder, major blood vessels and ureter,
- Cervix injury,
- Bleeding,
- Gas embolism,
- Pain,

- Intrauterine Adhesions,
- Infections,
- Vasovagal reactions.

Warnings

- The LiNA Versascope™ Outflow Cannula is supplied sterile using an ethylene oxide sterilization. Carefully inspect the packaging for any damage prior to use. Do not attempt to use the device if product sterile barrier or its packaging is damaged.
- Hysteroscopy procedures performed during pregnancy can result in miscarriage. Suspicion of pregnancy should suggest a pregnancy test prior to the performance of diagnostic hysteroscopy.
- When hysteroscopic instruments and accessories from different manufacturers are employed together in a procedure, verify compatibility prior to initiation of the procedure.
- Only fluid distension media can be used with the Schaft (and therefore with Outflow Cannula). LiNA Versascope™ Schaft is not applicable for hysteroscopy using CO₂ as distension medium.
- Strict fluid intake and output monitoring is required. Excessive intravasation of distention fluid can lead to fluid overload.
- Failure to follow all applicable instructions or any warnings or precautions may result in significant injury to the patient, physician or attendants and may have an adverse effect on the outcome of procedures performed. LiNA Versascope™ Outflow Cannula should not be used for any other purpose than intended function.
- Refer to appropriate electro-surgical system user manual indications and instructions to ensure that all safety precautions are followed.
- A thorough understanding of the principles and techniques involved in laser, electro-surgical, and ultrasonic procedures is essential to avoid shock and burn hazards to both patient and medical personnel and damage to the device and other medical instruments. Ensure that insulation or grounding is not compromised. Do not immerse electro-surgical instruments in liquid, unless the instruments are specifically designed and labeled to function in liquid.

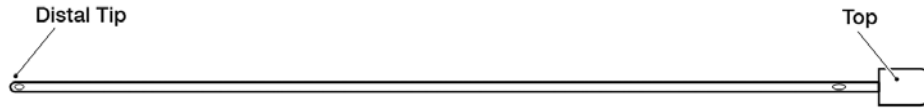- For single use only. Do not reuse, reprocess or re-sterilize the LiNA Versascope™ Outflow Cannula. Any reprocessing may impede the functions of the device. Reusing single use devices may also increase the risk of cross contamination which may lead to infection or transmission of bloodborne pathogens to patients and users. Attempts to clean the device results in risk of device malfunction.
- Do not use past expiration date.
- Do not use device if it has been exposed to non-sterile surfaces before procedures.
- Carefully inspect the device before use. If any damage or defect is identified, dispose of the device.
- No modification of this equipment is allowed.

Development Group for Management of Fluid Distension Media in Operative Hysteroscopy. BSGE/ESGE guideline on management of fluid distension media in operative hysteroscopy. Gynecol Surg. 2016; 13(4): 289–303),

- In case of elderly women with renal and cardiovascular problems, the upper fluid deficit levels should be 750 mL and 1500 mL, respectively. Therefore, careful monitoring of the fluid deficit during hysteroscopy is necessary (AAGL Practice Report: Practice Guidelines for the Management of Hysteroscopic Distending Media. J Minim Invasive Gynecol. 2013; 20(2): 137–148).

How supplied

The LiNA Versascope™ Outflow Cannula is supplied sterile in a Tyvek pouch for single use only. The device will remain sterile until the expiration date as long as the package is unopened and undamaged.

Device can be used in office, hospitals and ambulatory surgical suite procedure room by trained professionals.

Instructions for use

1. Visually inspect for breaches of sterile barrier system prior to use. Open the package and remove the LiNA Versascope™ Outflow Cannula using aseptic handling and transfer methods. Using sterile technique, remove the instrument from the package. To avoid damage, do not flip the instrument into the sterile field.
2. Insert the distal tip of the LiNA Versascope™ Outflow Cannula into the instrument port of the LiNA Versascope™ Sheath.
3. Fully insert the shaft of the LiNA Versascope™ Outflow Cannula into the LiNA Versascope™ Sheath until the top of the LiNA Versascope™ Outflow Cannula is flush with the Sheath housing.
4. When outflow is desired, open the roller clamp on the aspiration tube of the LiNA Versascope™ Sheath. Regulate desired distension and

outflow by adjusting the roller clamp on the outflow tube of the LiNA Versascope™ Sheath.

5. Hysteroscope AlphaScope should be reprocessed and disposed of according to instructions in the hysteroscope IFU. Sheath and Outflow Cannula should be disposed after use in accordance with hospital biohazard policy. Resterilisation the Sheath and Outflow Cannula devices may compromise the integrity of the device and lead to serious surgical consequences.

WARNING: Failure to open clamp on the outflow tube of the Sheath device and/or failure to use the Outflow Cannula may result in overdilation of the uterus or excessive intravasation of fluid.

Storage

No special storage conditions required. Do not use after expiration date.

Reporting

Any serious incident that has occurred in relation to the device should be reported to the LiNA Medical ApS and the competent regulatory authority of the country in which the user and/or patient is established.

CE marking 2022.

Country of origin: Poland.

Precautions

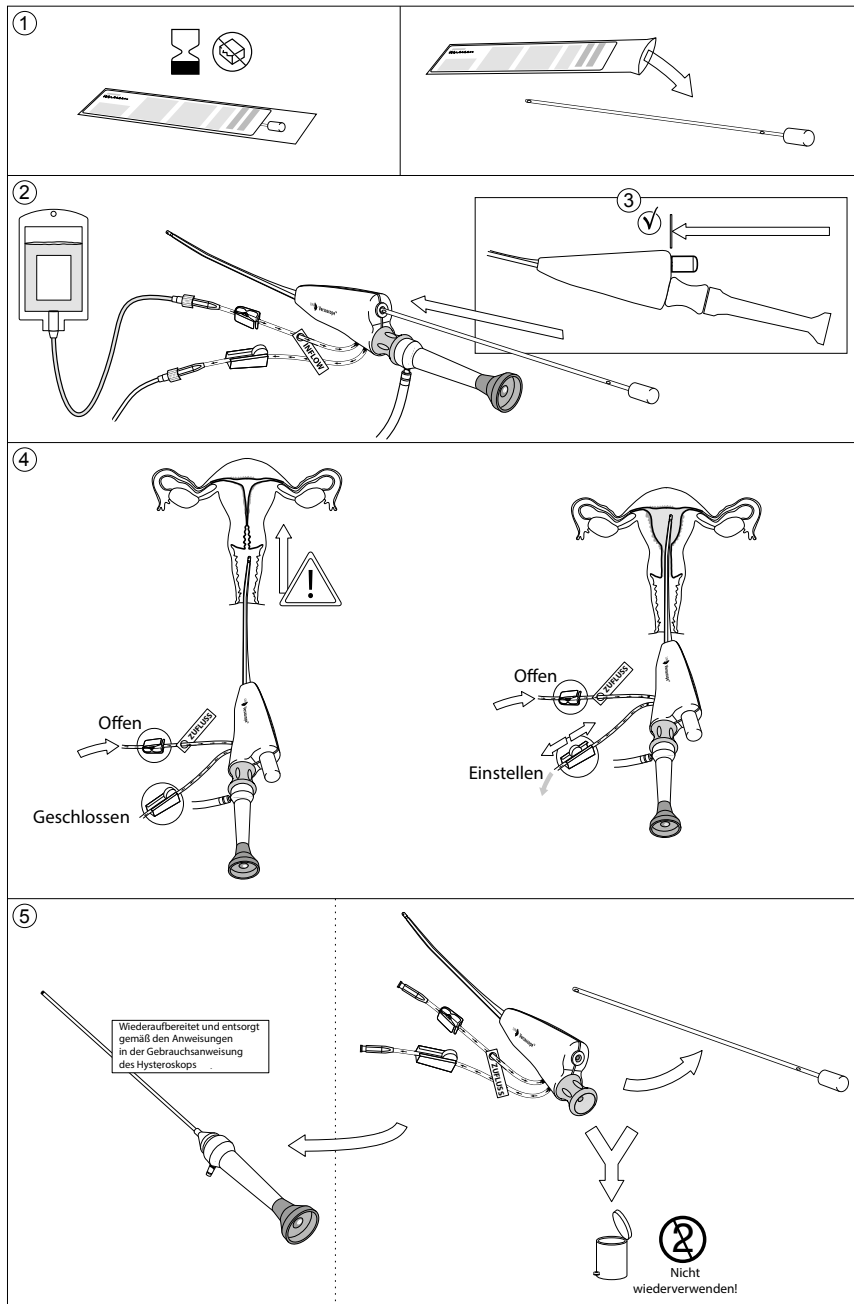
- After use, dispose of product and packing in accordance with hospital, administrative and/or local government policy.
- Only persons having adequate training and familiarity with hysteroscopy should perform hysteroscopic procedures. Consult the IFU of the compatible hysteroscope as well as medical literature relative to techniques, complications, and hazards prior to performance of any hysteroscopic procedure.
- Vaginal ultrasonography before hysteroscopy may identify clinical conditions that will alter patient management. Intrauterine distension can usually be accomplished with pressures in the range of 35–75 mmHg.
- Unless the systemic blood pressure is excessive, it is seldom necessary to use pressures greater than 75–80 mmHg.
- A fluid deficit (difference between the volume of fluid infused into the uterus and the volume of fluid evacuated from the uterus) should not exceed 1000 mL for hypotonic (glycine, mannitol, sorbitol) and 2500 mL for isotonic solutions (Umranikar S, Clark TJ, Saridogan E, et al. British Society for Gynaecological Endoscopy /European Society for Gynaecological Endoscopy Guideline

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